

GENICON®

Product Catalog 2013

New Ideas for Modern Healthcare





Suction Irrigation p. 26



Endoscopic Clip Applier p. 24



Open Clip Applier p. 22



Disposable Monopolar p. 20

GENICON Surgical Kits p. 18



Veress Needle / Specimen Retrieval Bag p. 16



Palm Bladeless p. 4

Bladeless Trocar p. 6

Shielded Trocar p. 8

Pyramidal Trocar p. 10

Hasson p. 12

Cannula p. 14



Palm Bladeless Trocar



SOCIETY OF
LAPAROENDOSCOPIC SURGEONS

- Ergonomic wing grip fits comfortably into your palm while its low profile provides comfort and stability
- Industry-first handle design, made from a single polycarbonate mold providing sturdiness and reliability with trocar placement
- Lean design reduces cost to produce product, with the savings passed on to the customer, a perfect option for budget conscious healthcare systems



- The distinctive shape of the bladeless tip separates and dissects tissue rather than cutting it, significantly reducing wound defect size, and minimizing insertion force



Palm Bladeless Trocar & No Cannula

Diameter	Length	Qty/Box	No Cannula
5mm	70mm	10	100-005-050
5mm	100mm	10	100-005-150
7/8mm	100mm	10	---
10mm	100mm	10	100-010-150
12mm	100mm	10	100-012-150

Palm Bladeless Trocar & One Cannula

Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula
5mm	70mm	5	100-005-051	100-105-051
5mm	100mm	5	100-005-151	100-105-151
7/8mm	100mm	5	---	100-108-151
10mm	100mm	5	100-010-151	100-110-151
12mm	100mm	5	100-012-151	100-112-151

Palm Bladeless Trocar & Two Cannula

Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula
5mm	70mm	5	100-005-052	100-105-052
5mm	100mm	5	100-005-152	100-105-152

Bladeless Visual Tip Trocar



SOCIETY OF
LAPAROENDOSCOPIC SURGEONS



- Universal handle can be used as a “Pistol Grip” when using a scope, or held as a traditional trocar handle. This grip versatility offers the surgeon a more controlled entry.



- High-visualization bladeless tip allows for unfettered visualization during insertion
- Injuries due to blind entry and a sharp blade are virtually eliminated with the use of the GENICON® Bladeless Tip Trocar System
- The distinctive shape of the bladeless tip separates and dissects tissue rather than cutting it, significantly reducing wound defect size

Bladeless Visual Tip Trocar


Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula	No Cannula
5mm	70mm	5	200-005-051	200-105-051	200-005-050
5mm	100mm	5	200-005-151	200-105-151	200-005-150
5mm	150mm	5	---	200-105-251	200-005-250
7/8mm	100mm	5	---	200-108-151	200-008-150
10mm	100mm	5	200-010-151	200-110-151	200-010-150
10mm	150mm	5	---	200-110-251	200-010-250
12mm	100mm	5	200-012-151	200-112-151	200-012-150
12mm	150mm	5	---	200-112-251	200-012-250
15mm	100mm	5	---	200-115-151	200-015-150
15mm	150mm	5	---	200-115-251	200-015-250

Bladeless Visual Tip Trocar & Two Cannula

Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula
5mm	70mm	5	200-005-052	200-105-052
5mm	100mm	5	200-005-152	200-105-152
5mm	150mm	5	---	200-105-252
7/8mm	100mm	5	---	200-108-152
10mm	100mm	5	200-010-152	200-110-152
10mm	150mm	5	---	200-110-252

Shielded Trocar

- The GENICON® Shielded Tip Trocar system incorporates a unique wedge style design, only the distal tip of the trocar provides a cutting action. The lateral sides of the blade are designed to dissect tissue rather than excise it. This results in a smaller defect size, and potential for reducing incidence of retroperitoneal herniation

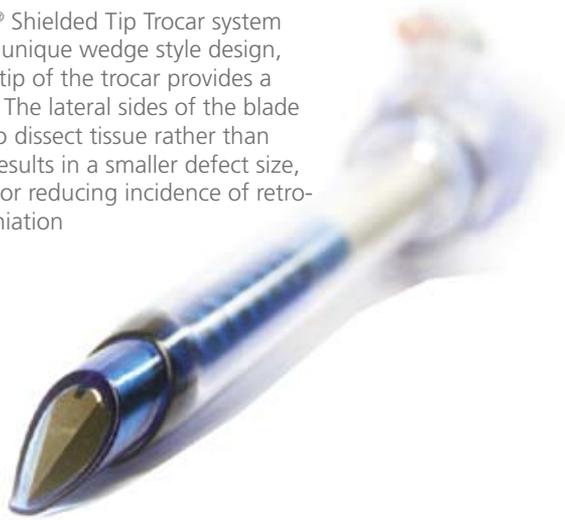

- Our Shielded Tip Trocar Systems are designed to be versatile, operationally flexible, and clinically adaptive to all surgical needs
- Ergonomic handle design allows for a firm, controlled grip, while minimizing hand fatigue



Shielded Tip Trocar & One Cannula

Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula
5mm	70mm	5	100-005-091	100-105-091
5mm	100mm	5	100-005-191	100-105-191
10mm	100mm	5	100-010-191	100-110-191
12mm	100mm	5	100-012-191	100-112-191

Shielded Tip Trocar & Two Cannula


Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula
5mm	70mm	5	100-005-092	100-105-092
5mm	100mm	5	100-005-192	100-105-192
10mm	100mm	5	100-010-192	100-110-192

Shielded Tip Trocar & No Cannula

Diameter	Length	Qty/Box	No Cannula
5mm	70mm	10	100-005-090
5mm	100mm	10	100-005-190
10mm	100mm	10	100-010-090
12mm	100mm	10	100-012-090

Pyramidal Trocar

- Industry-first handle design, made from a single polycarbonate mold, providing sturdiness and reliability with trocar placement
- Ergonomic handle design allows for a firm, controlled grip, while minimizing hand fatigue

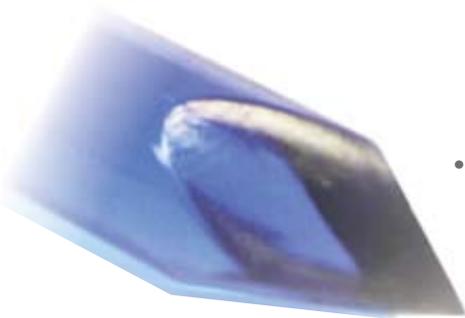

- Plastic-tipped Pyramidal Trocar uses the same cutting and insertion force as a comparable metal-tipped trocar

Pyramidal Trocar & One Cannula

Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula
5mm	70mm	5	100-005-001	100-105-001
5mm	100mm	5	100-005-101	100-105-101

Pyramidal Trocar & Two Cannula

Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula
5mm	70mm	5	100-005-002	100-105-002
5mm	100mm	5	100-005-102	100-105-102



Pyramidal Trocar & No Cannula

Diameter	Length	Qty/Box	No Cannula
5mm	70mm	10	100-005-000
5mm	100mm	10	100-005-100

Hasson



- The GENICON® Hasson offers an industry first ergonomic palm grip for more comfortable insertion

- Suture anchor wing assembly can be adjusted to the surgeon's desired depth using a single hand, while the unique suture anchor design allow for a quick and easy mode of securing the port in place
- Soft and malleable silicone sleeve conforms to skin and abdominal tissue, creating a better seal





Hasson

Diameter

12mm

Length

100mm


Qty/Box

5

Smooth Cannula

400-012-131

Cannula

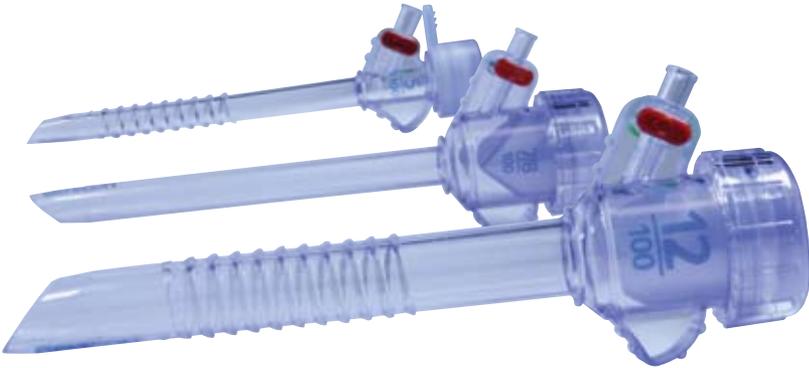

- GENICON Single-Use Cannula are available in a variety of configurations, and feature angled and beveled tips
 - Proprietary helical thread design provides a secure and reliable fixation to the abdominal wall
 - Intuitive push on / off button for CO₂ flow
 - Cannula cap is easily removed for specimen retrieval or rapid desufflation.



- Primary Valve will maintain seal with a variety of instrument sizes, from 4.7mm - 15mm.
- Unique geometry of the primary valve reduces instrument drag and resists tearing.

Cannula - No Trocar

Diameter	Length	Qty/Box	Threaded Cannula	Smooth Cannula
5mm	70mm	10	100-005-901	100-105-901
5mm	100mm	10	100-005-911	100-105-911
5mm	150mm	10	---	100-105-921
7/8mm	100mm	10	---	100-108-911
10mm	100mm	10	100-010-001	100-110-001
10mm	150mm	10	---	100-110-201
12mm	100mm	10	100-012-001	100-112-001
12mm	150mm	10	---	100-112-201
15mm	100mm	10	---	100-115-001
15mm	150mm	10	---	100-115-201

Veress Needle




- The GENICON® Veress Needle provides a high quality and cost effective means of insufflating the abdominal cavity
- Very low insertion force
- Intuitive push on / off button for CO₂ flow
- Tactile click upon entry into the abdominal cavity
- Ergonomic wings provide firm grip during insertion

Veress Needle

Description	Working Length	Qty/Box	Product Code
Veress Needle (Short)	120mm	25	130-120-001
Veress Needle (Long)	150mm	25	130-150-001

Specimen Bag

- Easy to use retrieval system designed to temporarily contain specimen and facilitate its removal from the patient during laparoscopic surgery
- Resilient, water-tight nylon bag material ensures confidence in specimen removal
- Flexible metal ring holds the bag open with 100mm opening for the Medium bag, 110mm opening for the Medium / Large bag, and 170mm for the Large bag
- Medium and Medium / Large bags may be easily introduced into a 10mm or larger cannula, Large bag may be introduced into a 12mm or larger cannula
- Cost effective and latex free
- Bag opens fully, allowing for easy specimen capture

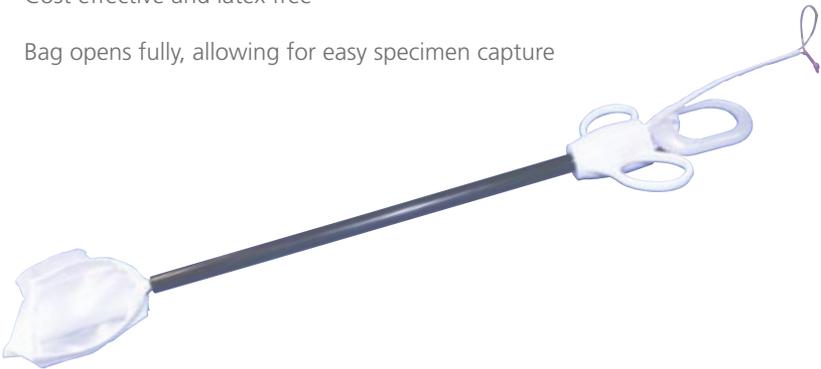

Specimen Retrieval Bag

Description	Bag Size	Qty/Box	Product Code
Medium Specimen Retrieval Bag	100mm x 150mm	10	550-000-003
Medium / Large Specimen Retrieval Bag	110mm x 200mm	10	550-000-004
Large Specimen Retrieval Bag	170mm x 240mm	10	550-000-005

GENICON Kits



- The GENICON Surgical Kit is available in several configurations, and contains the necessary array of trocars and insufflation needle to perform many common laparoscopic surgical procedures.
- New kits may be formed as demand for them is created; ask your territory sales representative for additional kit configurations that may be available.



GENICON Surgical Kit

GENICON, in an effort to meet the cost containment requirements of modern day healthcare, offers a vast array of Custom Procedure Kits. In partnership with Customers, GENICON will create Custom Procedure Kits, to meet the needs of any Laparoscopic Surgical Specialty, ranging from General, Gynecologic, Urologic, and Pediatric surgery.

Disposable Monopolar



Curved Metzenbaum Scissors
(Gray)



Straight Scissors
(Gray)



Fully Insulated Scissors
(Gray)



Curved Maryland Dissector
(Dark Green)



Atraumatic Fenestrated Grasper
(Blue)



Babcock Grasper
(Orange)



Maxi Grip Fenestrated Grasper
(Light Green)



- All GENICON® 33cm Disposable Monopolar Instruments feature color-coded rotation knobs to easily differentiate jaw configuration
- Instrument shaft rotation handle freely rotates the instrument tip 360° and is not limited by the jaw actuation



- Easy to use ratcheting handle function activated with a single-finger switchable toggle

Disposable Monopolar

Tip Description	Length	Qty/Box	Ratcheting	Non-Ratcheting
Curved Maryland Dissector	33cm	10	533-005-950	533-005-905
Atraumatic Fenestrated Grasper	33cm	10	533-005-920	533-005-915
Babcock Grasper	33cm	10	533-005-925	533-005-930
Maxi Grip Fenestrated Grasper	33cm	10	533-005-945	533-005-940
Curved Metzenbaum Scissors	33cm	10	---	533-005-910
Straight Scissors	33cm	10	---	533-005-935
Fully Insulated Scissors	33cm	10	---	533-105-910

Bariatric Disposable Monopolar

Tip Description	Length	Qty/Box	Ratcheting	Non-Ratcheting
Curved Metzenbaum Scissors	45cm	6	---	545-005-910

Open Clip Applier



- Distal-Tip closing enables control and precise placement.
- Smooth, tumbled tip reduces trauma to vessels
- Ligating Clips are sterile, single-patient-use clips produced from titanium



Straight Jaw



Jaw on angled applier positioned at 22°



Open Clip Applier

Jaw Description	Size	Qty/Box	Product Code	Clip Size
Straight Jaw	14.6cm	EA	330-005-000	Small
Straight Jaw	26.7cm	EA	330-005-001	Small
Straight Jaw	19.1cm	EA	330-005-002	Medium
Straight Jaw	26.7cm	EA	330-005-003	Medium
Angled Jaw	26.7cm	EA	330-005-004	M/L
Angled Jaw	26.7cm	EA	330-005-005	Large

6-Clip Cartridge

Size	Qty/Box	Product Code
Small	20	350-000-101
Medium	20	350-000-102
M/L	20	350-000-103
Large	20	350-000-104



Small



Medium




M/L



Large

Laparoscopic Clip Applier

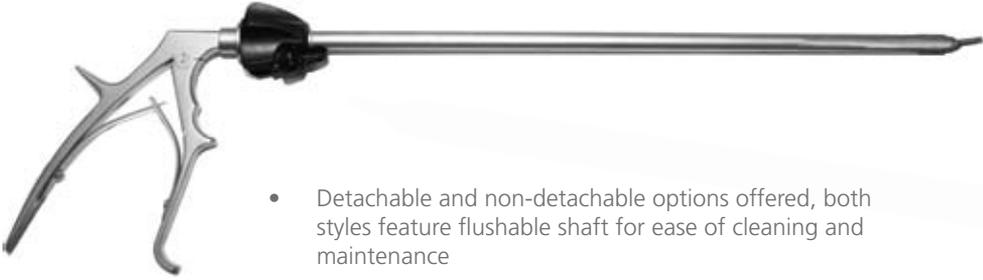

- Detachable and non-detachable options offered, both styles feature flushable shaft for ease of cleaning and maintenance



Single Action

- Available in single action and dual action tips



Dual Action



- Ligating Clips are sterile, single-patient-use clips produced from titanium

Clip Applier

Jaw Description	Size	Qty/Box	Product Code	Recommended Clip Size
Single-Action Jaw	5mm x 33cm	EA	350-005-000	Medium (White)
Dual-Action Jaw	10mm x 33cm	EA	350-010-000	M/L (Green)
Single-Action Jaw	10mm x 33cm	EA	350-010-001	M/L (Green)
Dual-Action Jaw (Detachable)	10mm x 33cm	EA	350-010-002	M/L (Green)
Single-Action Jaw (Detachable)	10mm x 33cm	EA	350-010-003	M/L (Green)

6-Clip Cartridge

Size	Qty/Box	Product Code
Small	20	350-000-101
Medium	20	350-000-102
M/L	20	350-000-103
Large	20	350-000-104



Small



Medium



M/L



Large


Suction Irrigation



- Probe tips are rounded, eliminating accidental lacerations to underlying tissue, and also provides an excellent tool for organ manipulation
- Removable rear port for easy access inside handpiece chamber

- Intuitively-colored buttons clearly distinguish between irrigation and suction functions

- Product is latex free and tubing is DEHP free

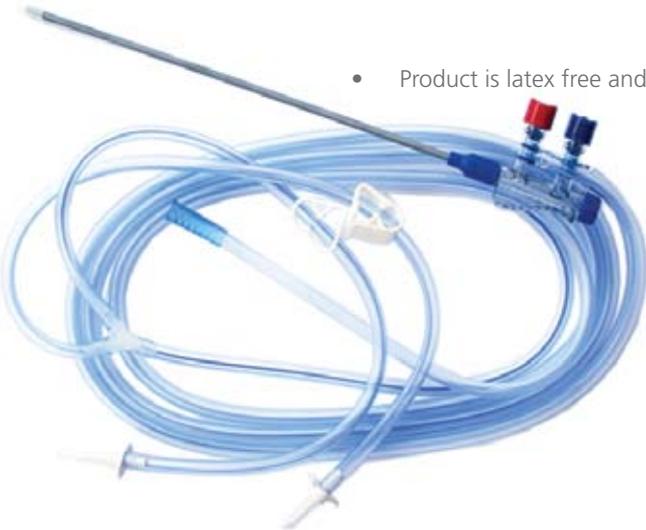

Disposable Suction Irrigation Tubesets

Description	Size	Qty/Box	Product Code
Gravity Fed Tubeset - No Probe (Single Clamp and Bag Spike)	---	10	700-000-005
Gravity Fed Tubeset	5mm x 28cm	10	700-005-005
Gravity Fed Tubeset	5mm x 33cm	10	700-005-010
Gravity Fed Tubeset (Single Clamp and Bag Spike)	5mm x 33cm	10	700-000-010
Gravity Fed Tubeset	10mm x 33cm	10	700-010-005
Gravity Fed Tubeset (2 Probes)	10mm x 33cm 5mm x 33cm	10	700-105-010
Powered Tubeset	5mm x 33cm	10	710-305-033
Powered Tubeset (Single Clamp and Bag Spike)	5mm x 33cm	10	710-300-000

Reusable Monopolar Probes

Description	Size	Qty/Box	Product Code
Spatula Probe	5mm x 33cm	1	707-005-133
J-Hook	5mm x 33cm	1	707-005-233
L-Hook	5mm x 33cm	1	707-005-333

Reusable Vent Probes

Description	Qty/Box	5mm x 33cm	5mm x 45cm	10mm x 33cm
Reusable Probe	1	707-005-033 (8-vent)	707-005-045 (8-vent)	707-010-033 (16-vent)

Disposable Vent Probes

Description	Qty/Box	5mm x 28cm	5mm x 33cm	5mm x 45cm	10mm x 33cm
8 Vent Probe	10	700-005-028	700-005-033	700-005-045	---
16 Vent Probe	10	---	---	---	700-010-033

Who we are ...

GENICON is recognized as an emerging leader in the design, production, and distribution of patented instrumentation focused exclusively on laparoscopic surgery. Our unique and innovative products meet the needs of surgeons and health care administrators while providing superior performance to the surgeon, a better outcome for the patient, and a lower cost to the healthcare provider.

Innovative Products

GENICON offers a full array of innovative products to meet your needs including: access ports, single use and reusable instrumentation and suction irrigation. GENICON products continue to lead the way in supporting the highest performance in laparoscopic technique. On multiple occasions, the Society of Laparoendoscopic Surgeons (SLS), has awarded the Innovation of the Year Award to GENICON for our Trocar Access Systems.



Olga Haberland

Director of International Business Development
Primary Phone + 1.407.657.4851 x304
Secondary Phone + 1.407.927.1743
USA Toll-Free 800.936.1020 x304
E-mail olh@geniconendo.com
Skype: [ohaberland](https://www.skype.com/user/ohaberland)
Languages Spoken: English, Spanish



Refaat a Moniem, MD

Middle East & Africa Market Manager
Primary Phone + 1.407.657.4851 x303
Secondary Phone + 1.407.413.0648
Fax +1.407.679.0880
E-mail rmoniem@geniconendo.com
Skype: [refaatmoniem](https://www.skype.com/user/refaatmoniem)
Languages Spoken: English, Arabic



Theresa Hoegstrom

U.S. Market Business Development Manager
Primary Phone + 1.321.305.8871
Secondary Phone +1.407.657.4851 x301
USA Toll-Free 800.936.1020 x301
E-mail tlw@geniconendo.com
Skype: [theresa.white.genicon](https://www.skype.com/user/theresa.white.genicon)



Amy Chen

Asia Business Development
亞洲市場部
Primary Phone + 1.407.657.4851 x302
Secondary Phone +1.407.741.3136
USA Toll-Free 800.936.1020 x302
E-mail achen@geniconendo.com
Skype: [amy.chen.genicon](https://www.skype.com/user/amy.chen.genicon)
亞洲市場部
Languages Spoken: English, Mandarin

For questions, comments, or orders, please email us at:
customerservice@geniconendo.com